

PORODO
LIFESTYLE



Porodo Lifestyle

Kids

Electric Toothbrush

SKU: PDLFSTK12



Table of Contents

Product Overview	2
Specifications	2
Safety Notes	3
Brush Head Installation	3
Schematic View	4
Product Features	5
Instructions for Use	5
Charging Guide	6
Care and Maintenance	6
FCC Statement	6
Disposal	7
Warranty	8
Contact Us	8

Product Overview

The Porodo Kids Electric Toothbrush is designed to make oral care fun and effective for children. It offers three customizable cleaning modes, allowing parents to select the best mode for their child's needs. With an IPX7 water-resistant rating, the toothbrush is safe to use in the shower or with water. The built-in smart timer ensures that children brush for the recommended duration, promoting better oral hygiene habits. On a full charge, the toothbrush lasts for up to 30 days, providing long-lasting use. Additionally, the toothbrush features Type-C fast charging, ensuring quick and efficient recharging. This toothbrush is perfect for children, offering both functionality and entertainment with its playful design.

Specifications

Input	5V
Rated Power	1.5W
Charging Interface	Type-C
Battery Capacity	350mAh
Charging Time	1H
Water Resistance Level	IPX7
Bristle Material	Dupont Bristles
Model Number	PDLFSTK12


Safety Notes

Some users may experience a tingling or pin-and-needle sensation when using the toothbrush for the first time. This is normal and will go away after a few uses.


- 1.** Do not open, repair, or alter the product without permission.
- 2.** Avoid knocking or dropping the product. If it is defective, stop using it.
- 3.** Ensure that the local voltage matches the marked voltage before connecting the power supply.
- 4.** Only use the toothbrush heads included in the box.
- 5.** If you have had oral or gum surgery in the past two months, consult your dentist before use.
- 6.** If your teeth bleed after using the toothbrush, or if bleeding continues for more than one week, consult your dentist.
- 7.** This toothbrush is for cleaning teeth, gums, and tongue only. Do not use it for other purposes.
- 8.** Children should be supervised to ensure proper use and should not play with the toothbrush.
- 9.** Do not place the toothbrush in high temperatures or expose it to direct sunlight.

Brush Head Installation

- 1.** Remove the Old Brush Head: Gently pull off the old brush head from the handle.
- 2.** Install the New Brush Head: Attach the new brush head onto the handle. A small gap of 0.5mm to 1mm may appear when installing the new brush head. This is normal and does not affect the performance of the toothbrush.


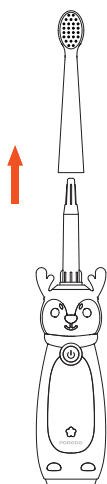


**1. Removing the
Old Brush Head**

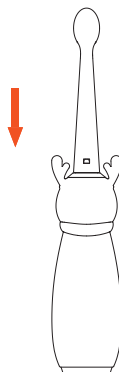

**2. Installing the
New Brush Head**

Schematic View

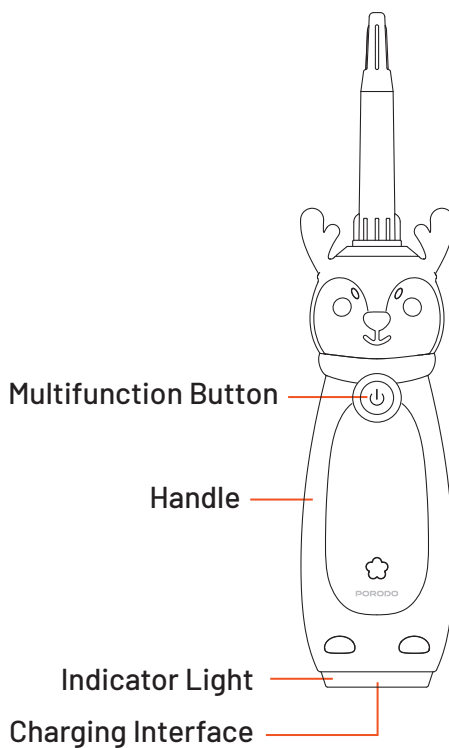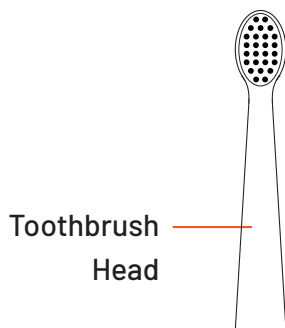

Product Features

- 1. Three Modes with Indicator Lights:**
 - a. Light Blue:** Sensitive mode for ages 3-6, designed for gentle brushing.
 - b. Blue:** Clean mode for ages 6-12, ideal for daily teeth cleaning.
 - c. Blue Flash:** Massager mode for ages 6-12, helps massage gums for improved oral health.
- 2. Waterproof Performance:** Rated IPX7, making the toothbrush safe for use in water, including the shower.
- 3. Long-Lasting Lithium Battery:** The battery lasts up to 30 days on a full charge, with a power shortage reminder and the convenience of USB charging for easy recharging.
- 4. Smart Mode Memory:** The toothbrush automatically stops after 2 minutes, features a 30-second timer reminder to move to different areas of the mouth, and encourages good brushing habits for effective oral care.

Instructions for Use

- 1.** Connect the USB cable to the power supply and charge the toothbrush.
- 2.** Attach the brush head securely. Wet the bristles and apply a small amount of toothpaste.
- 3.** Position the toothbrush with toothpaste against one section of your teeth.
- 4.** Gently move the brush head across each section of your teeth for 10 seconds, including the front, back, and chewing surfaces.
- 5.** To stop the device, press and hold the multifunction button for 3 seconds.
- 6.** Press the multifunction button to select one of the three modes: Sensitive, Clean, or Massage.
- 7.** Move to the next section when the 30-second auto reminder pauses. Repeat the steps for all four sections of your mouth. The toothbrush will automatically stop after 2 minutes.

Charging Guide

- 1.** If the red indicator flashes after use or the product cannot be turned on, it means the battery is low and needs charging.
- 2.** Connect the USB cable to a power supply to start charging. The indicator light will turn red while charging. Once fully charged, the indicator light will turn green.
- 3.** For first-time use or if the battery has not been charged for an extended period, it is normal for the indicator light to remain off for the first 5 minutes.

Care and Maintenance

a. Handle

- 1.** Clean the handle with a damp soft cloth and air dry it after removing the toothbrush head.
- 2.** Do not keep the handle in water for long periods.

b. Toothbrush Head

- 1.** Clean and air dry the toothbrush head after use.
- 2.** Do not press or pull the head hard, as it is easy to deform.
- 3.** To ensure optimal cleaning, replace the toothbrush head every 3 months.

FCC Statement

- 1.** This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:
 - (1)** This device may not cause harmful interference; and
 - (2)** This device must accept any interference received, including interference that may cause undesired operation.

2. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Reorient or relocate the receiving antenna.

Increase the separation between the equipment and receiver.

Connect the equipment to an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/TV technician for assistance.

RF warning statement:

The device has been evaluated to meet general RF exposure requirements. The device can be used in portable exposure conditions without restriction.

Disposal

This product must not be disposed of as unsorted household waste. It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations.

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at [**https://www.porodo.net/warranty**](https://www.porodo.net/warranty) and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
[**https://www.porodo.net/warranty**](https://www.porodo.net/warranty)

Contact Us

If you have any questions about this Privacy Policy, please contact us at:
[**info@porodo.net**](mailto:info@porodo.net)

Website: [**https://www.porodo.net/**](https://www.porodo.net/)

Service Support: [**support@porodo.net**](mailto:support@porodo.net)

Instagram: [**porodo**](#)